

CE
2460

EYECRYL PLUS CLEAR

Monofocal Hydrophilic Acrylic Aspheric
Foldable Intraocular Lens



biotech

VISIONCARE

EYECRYL PLUS CLEAR

Monofocal Hydrophilic Acrylic Aspheric
Foldable Intraocular Lens


FEATURES

- Aspheric optic surface
- 360° square edge
- Increased hydrophobicity because of surface modification hydrophilic IOL
- Less attractive to cell attachments & cell growth
- Surface modification helps in preventing the secondary cataract
- Exact refractive value of the IOL is printed along with diopter on the box
- 5° angulation ensures better contact with posterior capsule
- UV blocking Material

SPECIFICATIONS

MATERIAL	Hydrophilic Acrylic CQ UV		
OPTIC TYPE	Single Piece, 360° Square Edge with Aspheric Optic		
OPTIC SIZE	6.00 mm		
OVERALL SIZE	12.50 mm		
ANGULATION	5°		
ACD	4.96		
REFRACTIVE INDEX	1.462		
RECOMMENDED ULTRASOUND A-CONSTANT	SRK-T 118.00		
RECOMMENDED OPTICAL A-CONSTANTS	HAIGIS a0:1.030, a1:0.400, a2:0.100	HOFFER Q ACD: 5.13	Holl 1 : 1.37
	Holl 2 : 1.39	SRK - T 118.30	SRK - II 118.50
	Barrett: 1.52	Hill-RBF: 118.30	Olsen ACD: 4.16
DIOPTER RANGE	+5.0 D to +30.0 D (with 0.5 D steps)		
IMPLANTATION SITE	Capsular Bag		
STERILIZATION	Steam		
SHELF LIFE	4 years from date of manufacture		

Model: HSAS 600

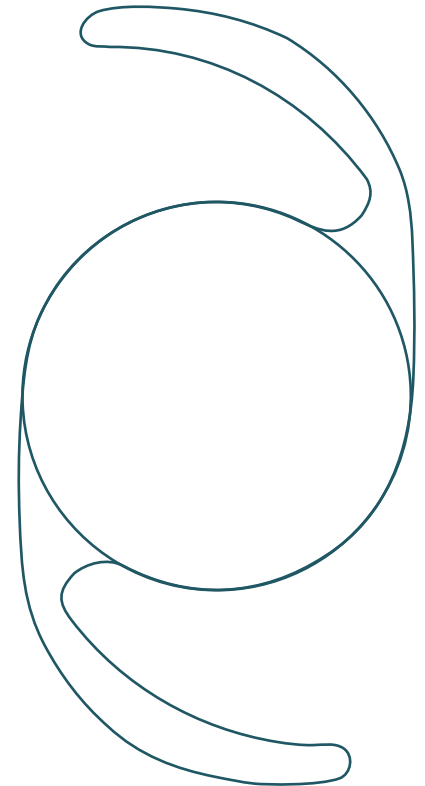

Scan For Biotech Connect App